Blister & Product Inspection



SPECTRA CL Colour Product Inspection



The colour product inspection system SPECTRA CL is used to detect cross-contamination and other errors such as size, position as well as broken products and overfilling. It is the ideal solution for common inspection tasks.

PRODUCT HIGHLIGHTS

Live data export 20 fps inspection speed



LDAP or OPC UA









scanware.de





DESCRIPTION

• Criteria: Colour, presence, size, shape, perimeter, position, broken products, overfilling, consecutive errors

• The system is fully compliant to 21 CFR Part 11

• scanware framegrabber for extra-fast evaluation


· Homogenous, reproducible and adjustable illumination - 16 levels available

Real-time operating system QNX®

• Soft- and hardware can be upgraded to SPECTRA HR if requirements change

· Also available in black & white for

simpler inspection tasks e.g. mono lines · Multitude of features depending on

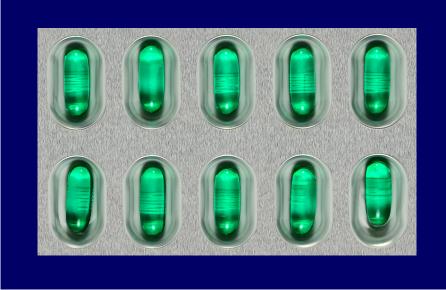

APPLICABLE OBJECTS

chosen software package



Soft gel

Hard gel

capsule



capsule

Oblong

Sugar-coated tablet

ON SITE


The software offers a number of features to optimise production on site: automatic teach-in and format test for quick start of production and adjusting to new batches. A number of features allow precise fine-tuning to the challenges of your product types. Finally, statistical features allow you to analyse product and the filling process.

scanware.de



GUI of software feature: Detection of upright objects

Using this feature enables you to minimise rejection of objects that are correct, but sat upright, without increasing tolerances

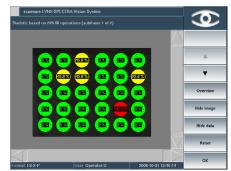

GUI of software feature: Pocket statistics This feature allows you to assess your filling for consecutive errors and carry the gained knowledge over to your filling machinery

SCANWARE HIGHLIGHTS



info@scanware.de

Any questions?

+49 6257 9352 0